



Pill Buddie

MEDICATION PLANNING GUIDE

Carousel loading reference card

Before You Begin

This quick reference shows how to load your PillBuddie carousel based on how many times per day you take medication.

Never turn the carousel by hand.

Doing so will void the warranty and damage the unit.

Starting a New Cycle

- Always start a new cycle with the covered HOME slot positioned over the tray.
- Perform the HOME cycle first if the carousel is not already aligned.

Loading Your Tray


- Follow the color-coded guide on the following panels to fill your tray based on how many times per day medication is needed.
- Always start loading at Slot 1 and work sequentially around the wheel.
- Confirm your exact daily dose count and times with your pharmacist or physician before setting up your schedule.

What's Inside

- 1 Dose/Day — 28-Day Supply
- 2 Doses/Day — 14-Day Supply
- 3 Doses/Day — 9-Day Supply
- 4 Doses/Day — 7-Day Supply

1 Dose Per Day


28-Day Supply

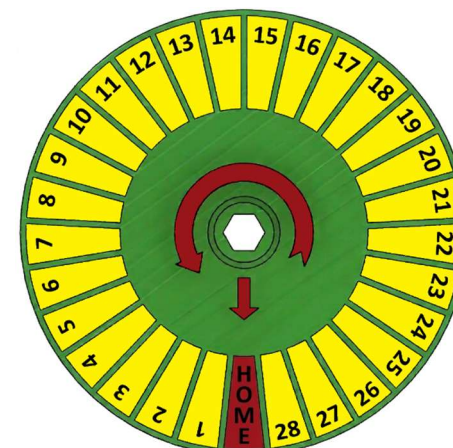

- Fill all 28 slots with your scheduled daily medication.
- Provides a full 28 days of dosing before a refill is needed.

2 Doses Per Day

14-Day Supply

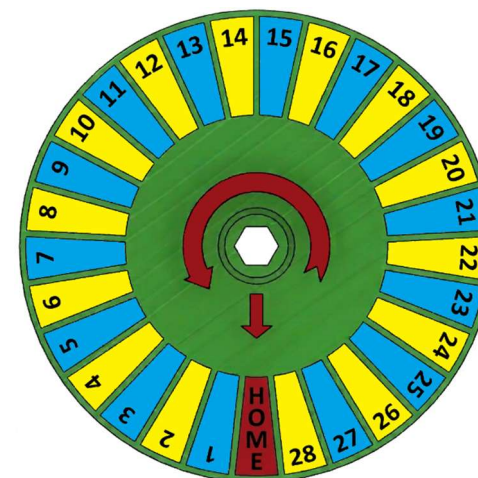

- Alternate medications in every other slot.
- **Yellow Slots** — Dose 1 (e.g., Morning)
- **Blue Slots** — Dose 2 (e.g., Evening)
- Provides 14 full days (2 weeks) of dispensing.

3 Doses Per Day

9-Day Supply



- Rotate medications across every 3 slots.
- **Red Slots** — Dose 1 (e.g., Morning)
- **Yellow Slots** — Dose 2 (e.g., Afternoon)
- **Purple Slots** — Dose 3 (e.g., Evening)
- Provides 9 full days of dosing, with 1 extra slot remaining.

4 Doses Per Day

7-Day Supply



- Rotate medications across every 4 slots.
- **Purple Slots** — Dose 1 (e.g., Morning)
- **Yellow Slots** — Dose 2 (e.g., Noon)
- **Red Slots** — Dose 3 (e.g., Evening)
- **Blue Slots** — Dose 4 (e.g., Bedtime)
- Provides 7 full days (1 week) of dispensing.

Need Help?

Not sure which schedule fits your routine? Reach out — we're happy to help.

Cleaning & Maintenance

- Wipe the exterior and tray with a soft, damp cloth.
- Never immerse in water or use harsh chemicals.

info@PillBuddie.ca

(902) 989-1807

pillbuddie.ca

Scan for support:



PILLBUDDIE

Healthy Tech

© PillBuddie Healthy Tech